


# Getein 1160 Immunofluorescence Quantitative Analyzer Installation Guide

- Check if instrument packing case is damaged before unpacking. If it is, please contact the carrier and our after-sales support.
  - Please well-preserve the packing case for future use.
  - The analyzer should be placed in a dry and clean area, which is free from dust, vibration, strong electromagnetic interference, direct sunlight, and large temperature difference.
1. Open the packing case, take out the analyzer, and loosen the thumb screw counterclockwise at the bottom of the analyzer.

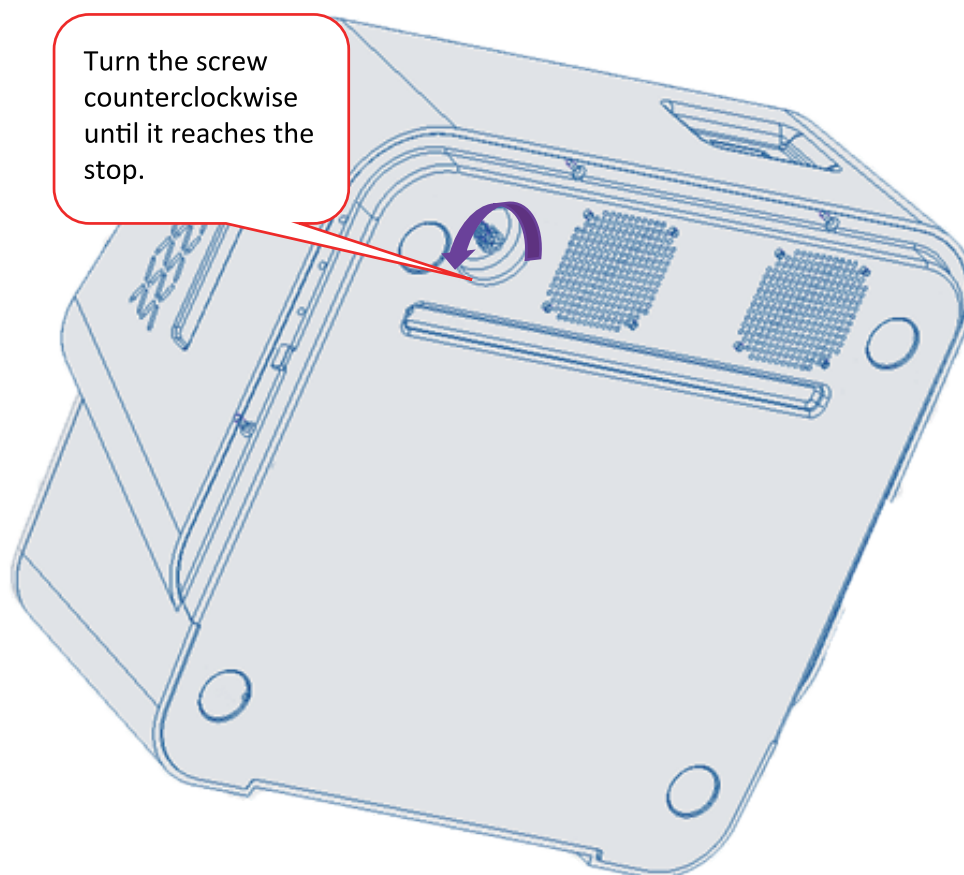

2. Put the analyzer on a flat platform, load printing paper into built-in printer, and connect the power plug with an electric outlet.
3. Please read the user manual for more details.