



## 1. Appearance Description

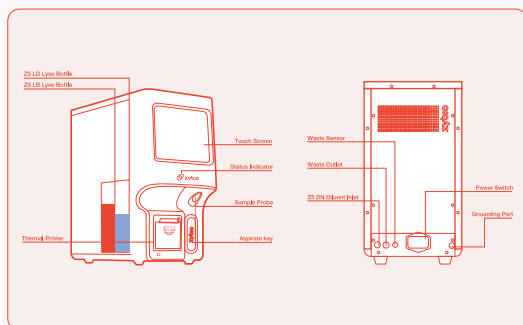

## 2. Daily Operating Process

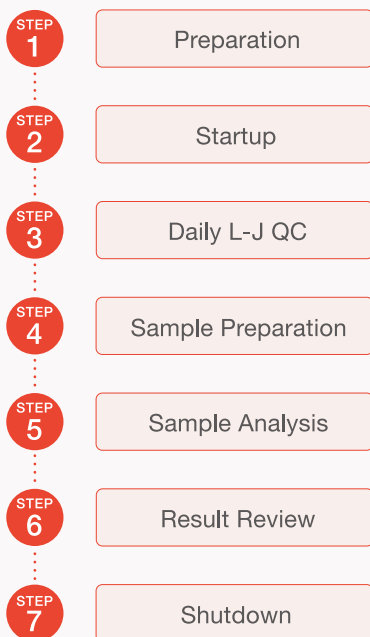


## 2.1 Preparation

- Check the waste bucket and make sure it's not full.
- Check the reagents to see if the reagents are expired or frozen.
- Check the thermal printer and make sure enough printer paper is installed.

## 2.2 Startup

- Press power switch at the back of the analyzer.
- Input **username** and **password**. (Default: Username: admin, Password: admin)

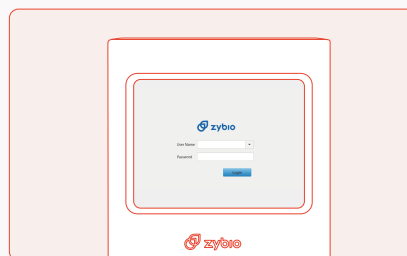

## 2.3 Daily L-J QC

- Tap the menu option "QC" => "L-J QC" => "Setting".
- Enter the "L-J QC setting" interface. Tap "New", or select a QC file without QC results.
- Enter the expiration date of the lot.
- Select "control level", "control type" and "QC mode" and set "QC sample ID".
- Enter the target and limits in the edit boxes according to the package insert of the lot of controls.

- Take out one level of control from the refrigerator and handle it as the control insert.
- Run control under the "QC" interface and review the "QC graph" or "QC table".

## 2.4 Sample Preparation

- The analyzer can run 3 types of samples: venous whole blood samples, capillary whole blood samples and prediluted samples.
- Venous whole blood: collect at least 0.5 mL EDTA-K<sub>3</sub> or EDTA-K<sub>2</sub> whole blood and run it within 15 min to 8 hours after being collected to ensure the accuracy of the results.
- Capillary whole blood: collect at least 100 µL whole blood and run it within 5 min to 2 hours after being collected to ensure the accuracy of the results.
- Blood samples shall be mixed as shown in the following figures.

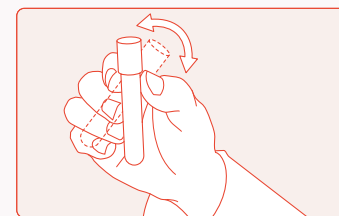


Figure 1 Upside down and mix at least 8 times

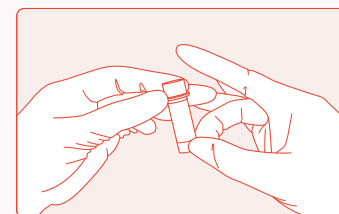


Figure 2 Flip and mix well (Used for capillary whole blood)

## 2.5 Sample Analysis

### 2.5.1 Venous Whole Blood Analysis



Tap “**Analysis**” and tap “**next sample**”.



Choose “**VWB-CBC**” or “**VWB-CBC**” or “**VWB-CBC+ DIFF**” analysis mode (default mode: “**VWB-CBC+DIFF**” )



Input patient information.  
Mix the sample as above 2.4.



Put the sample under the sample probe and press aspirate key.

### 2.5.2 Capillary Whole Blood Analysis



Tap “**Analysis**” and tap “**next sample**”.



Choose “**CWB-CBC**” or “**CWB-CBC+DIFF**” analysis mode.



Input patient information.  
Mix the sample as above 2.4.



Put the sample under the sample probe and press aspirate key.

### 2.5.3 Prediluted Blood Analysis



Tap “**Diluent**” and use a clean centrifugal tube to aspirate 480  $\mu$ L diluent.



Collect 20  $\mu$ L of capillary blood with a tube contained diluent. Mix the sample as above 2.4.



Run it within 5 min to 30 min after being collected to ensure the accuracy of the results.



Choose “**PD-CBC**” or “**PD-CBC+ DIFF**” analysis mode and input patient information.



Put the mixed sample under the sample probe and press aspirate key.

## 2.6 Result Review



Select “**Review**”.



Select relevant results to manage: “**Graph Review**”/ “**Trend graph**”/ “**Print**”/ “**Export**”/ “**Validate**”.



Tap and enter “**Graph Review**” Interface to see the parameter, histogram and scattergram.

## 2.7 Shutdown



Tap “**Shutdown**”.



Put Probe cleanser under the sample probe and press aspirate key.



Press the power switch after cleaning.

## 3. Replace Reagent

When LD/LB Lyse or DN Diluent is ran out.



Please tap “**Reagent**” => “**Replace Reagent**”.



Take out a new reagent and get the new RFID card. Put it on the RFID reader that can be automatically identified to get the specification of reagent.



Connect the new reagent.



Tap “**Apply**” and prime the reagent automatically.

## 4. Routine Maintenance

• When analyser shows any alarm, please tap “**Alarm**” information interface and tap “**Clear**”.If alarm still exists, please follow steps below.

• **[Background Abn./ Sample result is not good]** Tap “**Management**”=> “**Service**” => “**Maintenance**”=> “**Whole Device**”, then start probe cleanser soaking.

• **[HGB Blank Voltage Abn.]** Tap “**Management**”=> “**Setup**”=> “**System**” => “**Gain**” => “**4.5 V**”=> “**OK**”, then retest one sample.

• **[LB Lyse/LD Lyse/DN Diluent Empty]** Tap “**Management**” => “**Service**” => “**Maintenance**” => “**DN Diluent**” or “**LD Lyse**” or “**LB Lyse**”, then start to prime reagents.

• When the analyser reminds you to do the “**cleanser soak**”, perform it as instruction.

• If the analyser is not to be used over 1 week or needs transportation, perform “**Drain All**” as instruction.