



How to cancel “Reagent Expired” alarm

Sometimes when reagent is out of the expired date, the instrument will ask for replacing reagent while you want to test a sample. If the reagent performance is good enough and the user still want to continue to use it, you can follow the below steps:

1. Enter **Management** => **Setup** => **Reagent** menu.
2. Select “One key elimination” to activate the function and click **OK** to save.

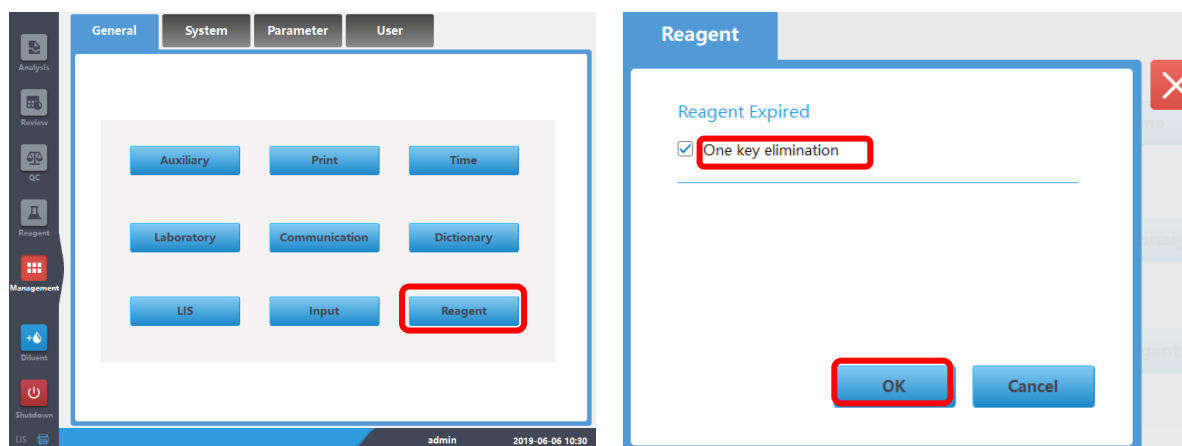

3. Then we can click the alarm area, on the pop-up dialog have one key clear alarm function. Click **clear** to cancel the alarm and close the dialog box, then we can continue use the reagent.

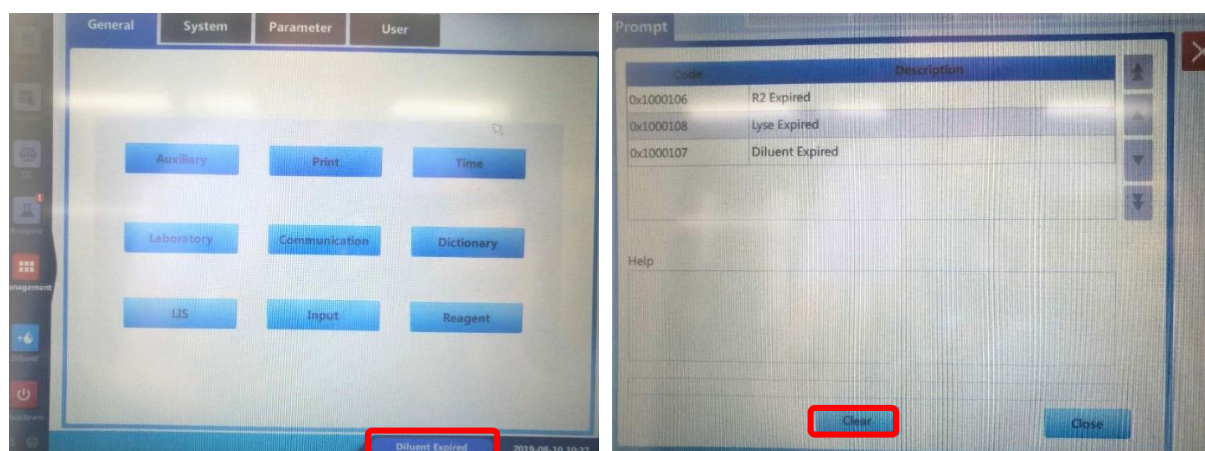

Alarm Area